WiseFido Radar

Quick start · mount it, then set it up in the owlCare app


Full guide & help
wisefido.com/install

In the box


① Radar · ② Ceiling bracket · ③ Wall bracket · ④ Power adapter · ⑤ Power cable (USB-C, 6 m) · ⑥ Cable clips · ⑦ Adhesive bases · ⑧ Steel pins — ⑦+⑧ mount the bracket without drilling.

1 Mount the radar

Ceiling · 2.3–3.0 m (7.5–10 ft) — middle of the room, above the bed if you can. Indicator light points away from the door's wall.

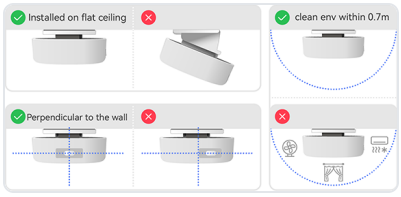

Bedroom wall · about 160 cm (5 ft 3 in) — at the head of the bed, facing into the room. Not on the wall the door is on.

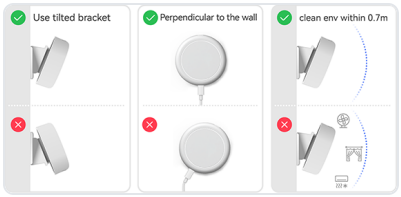


Bathroom · on top of the mirror (≈205 cm) or on the door frame (≈210 cm), facing into the room — **enter 210 cm**.

Height = floor to the surface the radar sits on. Keep fans, A/C and curtains 1 m (3 ft) away, heat lamps 50 cm (20 in).

2 Add it in owlCare

Get the app:

iPhone Android

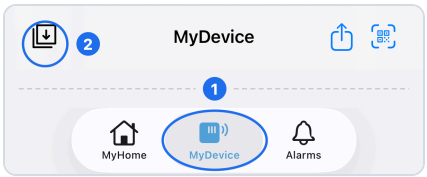

① **MyDevice** → ② **Add device** (top left) → **Scan QR code** (on the back of the radar) or **Nearby (Bluetooth)** → pick your radar → enter WiFi name and password → **Connect** → **Choose location**: the room it is in.

2.4 GHz WiFi only. Keep your phone next to the radar.

3 Set the install mode


MyHome → **Build** → tap your radar. On **1 Radar**:

- **Mode**: Wall or Ceiling — how you mounted it.
- **Height**: from step 1 (bathroom: 210 cm). Green ✓ = recommended.
- A wall longer than 5 m? Tick **Large Room**.
- Tap **Save**.

The top bar shows ✓ Wall or ✓ Ceiling. ⚠ or ✕ means the radar is not seated in its bracket — push it back in.

4 Mark the bed

Open **3 Bed** → **Add**. The person on the map is you, live.

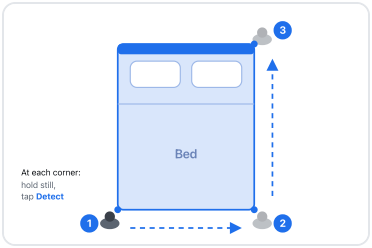

1. Walk to a corner of the bed and hold still, tap **Detect**.
2. Same for the next 2 corners, then tap **Save**.

5 Door and sofa Optional

- **Door**: learned automatically. To set it now: **2 Entrance** → **Add**, stand at each side of the doorway.
- **Long sofa**: **3 Bed** → **Long Sofa**, sit on each end.

6 Check it works

- Lie on the bed: in about a minute, owlCare shows heart rate and breathing.
- After a fall, alerts come after 60 s by default (adjustable).
- Restart WiFi setup: hold the radar button 10 s. More help: **wisefido.com/install**.

Not a medical device. It does not diagnose, treat or prevent any disease, and does not replace emergency services or personal care.

Warranty: 1 year from the date of purchase. Contact sales@wisefido.com.

⚠ **WARNING:** Cancer risk from exposure to lead. See www.P65Warnings.ca.gov.

Distributed by WiseFido LLC · 5171 Fox Meadow Dr, Highlands Ranch, CO 80130, USA · +1 (720) 694-0045 · Made in China

FCC ID: 2BD97HC2 · AI Care Sensor · Grantee: TsingLan Technology (Shenzhen) Co., Ltd

RF exposure: This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. It should be installed and operated with a minimum distance of 20 cm (8 in) between the radiator and your body.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: reorient or relocate the receiving antenna; increase the separation between the equipment and receiver; connect the equipment into an outlet on a circuit different from that to which the receiver is connected; consult the dealer or an experienced radio/TV technician for help.