
WiseFido Sleepad

Quick start · under the mattress, then set it up in the owlCare app

Full guide & help
wisefido.com/install

1 Position on the bed

Place the pad **50–60 cm (20–24 in) from the headboard**, centered across the bed, with at least **25 cm (10 in)** free on both sides.

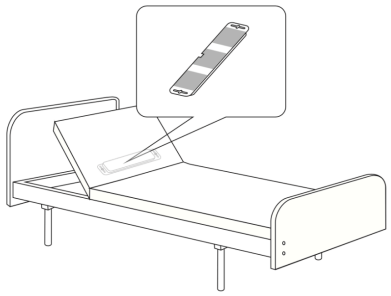


2 Under the mattress


Lay the pad flat on the bed base, **under the mattress**. Stick its back to the base with the 3M tape so it does not slide.

Not for water beds or air mattresses.

3 Connect the cables


Pad cable into the **Sensor Pad** port, power adapter into the **Power** port, then plug the adapter into a wall socket.

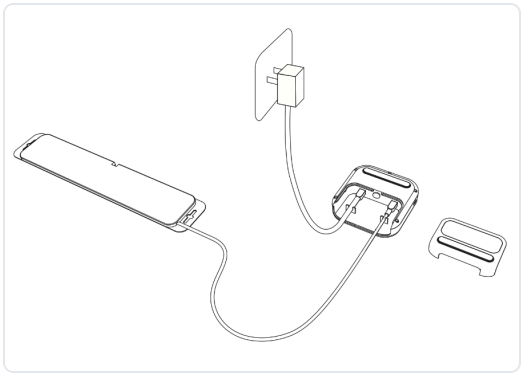

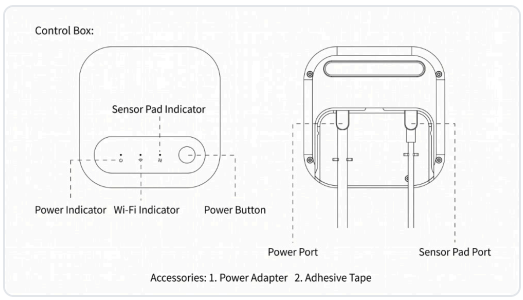

4 Add it in owlCare

Get the app: iPhone Android

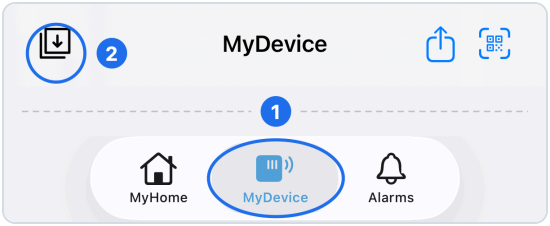

① **MyDevice** → ② **Add device** (top left) → **Scan QR code** (on the back of the control box) or **Nearby (Bluetooth)** → pick your Sleepad → enter WiFi name and password → **Connect** → **Choose location**: the bed it is on.

2.4 GHz WiFi only. Keep your phone next to the control box.

5 Check the lights

- All lights off or not red: working.
- **WiFi light red**: not connected to WiFi. Repeat step 4.
- **WiFi light flashing red**: WiFi is fine but the server is unreachable. Check your internet.
- **Sensor Pad (Link) light red**: pad not connected. Re-seat the pad cable.

FCC ID: 2ADIOBM8701-2 · Model: BM8701-2 · CMIIT ID: 2018DP2463

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: reorient or relocate the receiving antenna; increase the separation between the equipment and receiver; connect the equipment into an outlet on a circuit different from that to which the receiver is connected; consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Not a medical device. It does not diagnose, treat or prevent any disease, and does not replace emergency services or personal care.

Warranty: 1 year from the date of purchase. Contact sales@wisefido.com.

⚠ **WARNING:** Cancer risk from exposure to lead. See www.P65Warnings.ca.gov.

Distributed by WiseFido LLC · 5171 Fox Meadow Dr, Highlands Ranch, CO 80130, USA · +1 (720) 694-0045 · Made in China